



FlexiTab® S Manual tablet press

2075 0000 product line

Manual tablet presses for galenics and research.

With the FlexiTab® S you can now easily press a tablet with just a few milligrams and without a long set-up time. Pressing takes place manually or if desired electrically.

The simple, but well thought-through concept of the FlexiTab® S ensures you have plenty of flexibility in handling your granulates.

You can press both round and shaped tablets. You are optimally equipped for nearly all requirements with a pressure of up to 50 kN (5 tons) and a maximum tablet diameter of 20 mm.

The product-contacting parts are either made from stainless steel or pharmaceutical grade plastic. Thus the press can also be used for pharmacy production or for clinical samples.

Stamps of type EURO/TSM-B can be ordered directly from us. Optionally round, oblong, oval or free-form; with bevel, dividing groove or logo. Please get in touch with us.

Subject to changes. 03 -15



FlexiTab® S Manual tablet press


2075 0000 product line

Advantages of the FlexiTab® S:

- Quick and easy setting of the filling depth
- Quicker and easier format change without removal of the filling shoe
- Easy cleaning of the product-contacting parts
- Product-contacting parts of stainless steel or POM
- Pressing of very small quantities by "hand filling"
- Readily accessible press chamber
- Easy to transport



Pressurisation by manual hydraulic pump

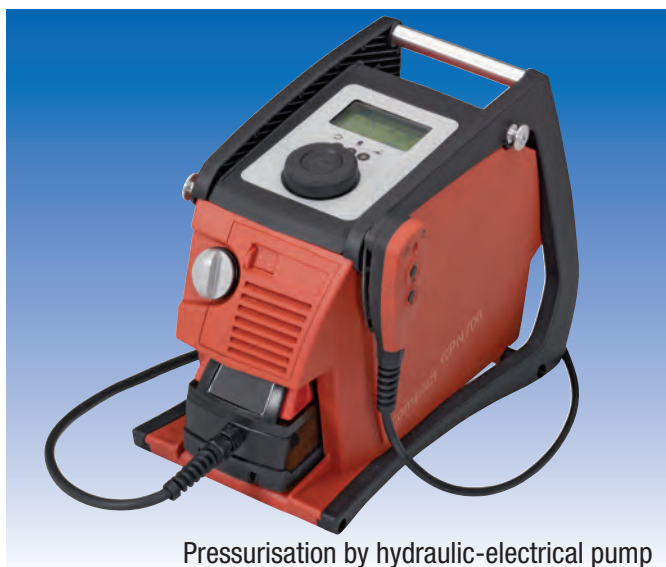

Pressurisation by hydraulic-electrical pump

Technical data:

Stamping force:	50 kN
Max. tablet diameter:	20 mm
Max. filling depth:	20 mm
Max. immersion depth, top stamp:	12 mm
Dimensions (W x H x D):	250 x 500 x 200 mm
Product-contacting parts:	V4A or Delrin®
Weight:	approx. 55 kg

Subject to changes. 03 -15